


MR BREAST SOLUTIONS MADE BY NORAS

Breast BI 7 MR Coil Set 1.5T & 3T

Exclusively for the new Siemens MRI scanner generation.

Patented lighting system
integrated

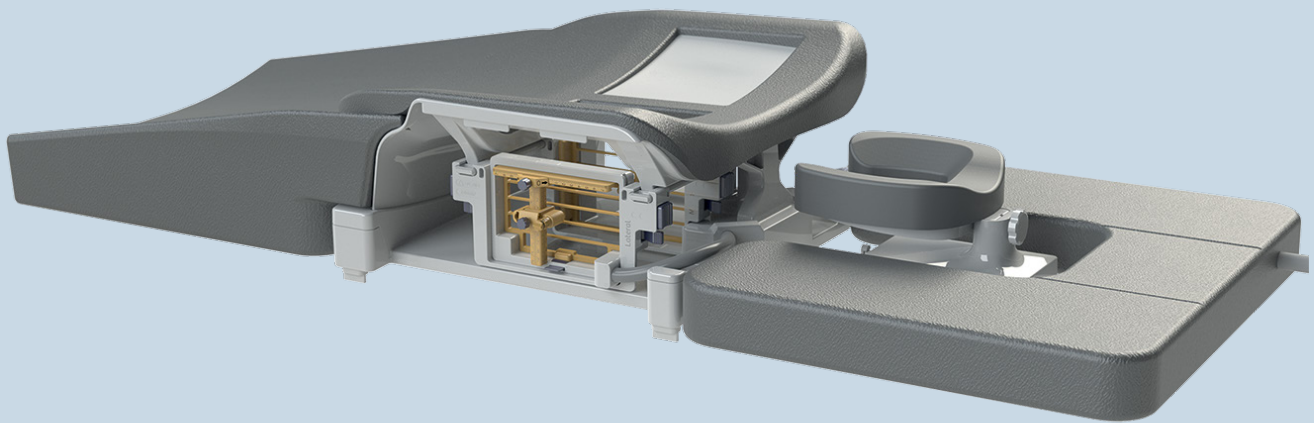


The **NORAS Breast BI 7 MR Coil Set** is an innovative solution designed to offer flexible use to meet the needs of the procedure to be performed in your clinic. Using the same set of comfort cushions, it can be combined with the Breast 18 imaging coil from Siemens. The rigid coil frame forms a high element-density, 7-channel „phased array” configuration for high resolution imaging. One sided medial and up to two-sided lateral biopsy access is possible.

Advantages of the Breast BI 7 MR Coil Set:


- Breast biopsy and diagnosis solution, also for very large breast volumes (11 liters)
- Simple and wide access for breast biopsy (lateral and medial)
- Excellent image quality with high homogeneity
- Very good imaging of the axilla
- The integrated lighting system (LED) with a contactless proximity sensor supports the biopsy workflow
- Medio-lateral or cranio-caudal fixation of the breasts possible for diagnostic imaging

The Post & Pillar biopsy system is designed for targeted needle guidance, which can be moved both left-right and anterior-posterior. Additionally, it supports angulated biopsies. As an alternative, disposable biopsy grids* can be used to perform biopsies also at numerous lesions at the same time. Disposable adapters* for all compatible biopsy needles and vacuum guns are available directly from NORAS.

* distributed by NORAS

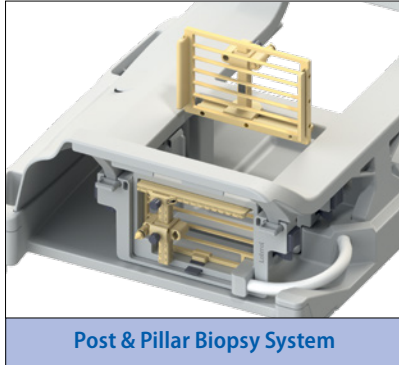

Post & Pillar Biopsy System

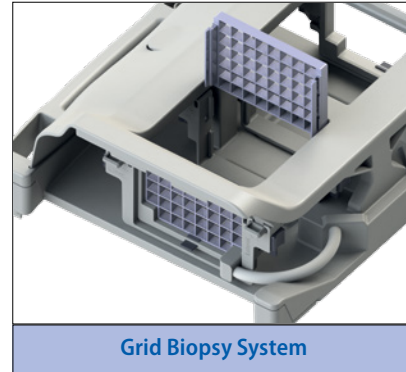


Grid Biopsy System

Specially developed fixation units for medio-lateral or also cranio-caudal fixation of the breast are available as an option to avoid motion artifacts and reduce scan time.



Cranio-caudal Fixation

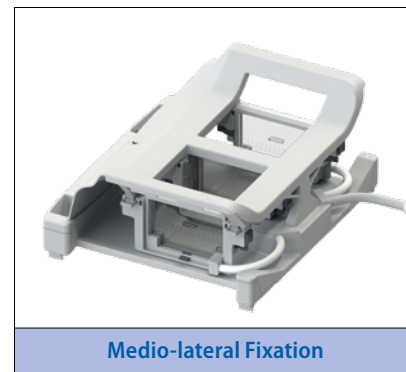

Medio-lateral Fixation

Scanner-Compatibility

If you have any questions regarding the compatibility of the Breast BI 7 MR Coil Set 1.5T & 3T, please kindly contact your Siemens sales representative.

NORAS – Your partner for innovative solutions in MR Imaging

Since more than 30 years NORAS MRI products GmbH is a reliable partner for high resolution imaging in MR-based diagnostics and intervention. Kindly get in touch with us in order to learn more about the application areas of our solutions for your MR-based tasks.

NORAS MRI products GmbH
GERMANY
 Leibnizstrasse 4
 97204 Hoechberg
 Tel.: +49 (0) 931 2 99 27 0
 Fax: +49 (0) 931 2 99 27 20
 info@noras.de
 www.noras.de



NORAS MRI products GmbH
USA
 Mr. John Wilkie
 Tel.: +1 (414) 350 4433
 sales-usa@mri-products.com
 www.mri-products.com